

PAJUNK®



NRFit®

Products & Services for *Regional Anaesthesia* Connectors
acc to ISO 80369-6



MADE IN GERMANY 


NRFit® – The New Standard for Regional Anaesthesia

In order to avoid the risk of incorrect connections and injections in the future, the ISO 80369-series of standards was issued in 2016. The new ISO norm replaces the overarching LUER standard and defines a separate, incompatible connection for each of the five application areas within the standard series: 80369-2 for Breathing Systems and Driving Gases, 80369-3 (**ENFit®**) for Enteral and Gastric, 80369-5 for Limb Cuff Inflation, 80369-6 (**NRFit®**) for Neuraxial procedures and 80369-7 (LUER) for Intravenous applications. As the leading manufacturer in Regional Anaesthesia, PAJUNK® offers its entire product portfolio for Neuraxial, major Regional Anaesthesia, Wound Infiltration as well as lumbar puncture applications in accordance with ISO 80369-6. This brochure shows you the product categories available in the new standard and our service package for a successful conversion to NRFit.


NRFit at a glance – what changes?



LUER → NRFit®

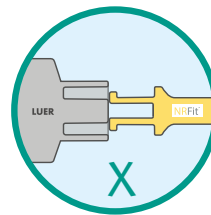

NRFit® → LUER

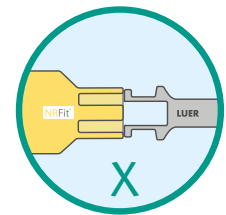

New connector


All connectors of the needles and components are affected by the changes

- NRFit connectors are **20% smaller and 3 mm longer** than Luer connectors. The NRFit slip connector is equipped with an additional outer cone.

Avoidance of misconnections

The probability of misconnection or misinjection is minimised by the new series of standards

- NRFit connectors can't be connected to LUER connectors
- For a better distinguishability all NRFit connectors are colour-coded in yellow

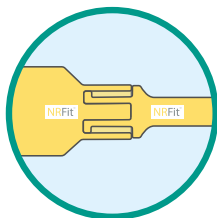

Limited compatibility

The compatibility of the NRFit connectors is limited to the neuraxial and major Regional Anaesthesia applications defined under ISO 80369-6

Change in product name

Products that meet this standard (ISO 80369-6) are given the additional appendix "NRFit"

Fit for NRFit

The PAJUNK® All-Round Service for a Successful Conversion



Ready to use: Standard and Customised NRFit Sets

To make the LUER to NRFit transition as easy as possible, Pajunk offers complete standard sets with additional NRFit components (syringes, draw-up needles, straws, etc.) as well as customised set options which allow for full customisation and flexibility. No matter which NRFit set option you prefer, we will find the right configuration for your individual application to support your needs and the existing facility processes.



Identification: Clear Labelling & Product Identification

The yellow NRFit markings on products & packaging allow for easy identification, while reducing confusion.

Counseling: Planning, Partnering & Execution

Various areas such as procurement, storage and training, maintenance and application are affected by the conversion process.

Pajunk supports you with the necessary information, helps you to correctly assess the risks of the conversion process and accompanies you in the development of the conversion concept with various services around the following aspects:



Procurement management

- Change of article numbers
- Needs analyses for the respective application areas

Storage

- Separate storage of the individual connector types
- Demand-related storage of all required parts

Medical application

- Application training with the new connectors
- Determination of all affected processes under Regional Anesthesia

The Pajunk NRFit product range:

Spinal/Epidural/CSE, Nerve blocks, Lumbar puncture, Wound infiltration

NEURAXIAL ANAESTHESIA

Spinal Anaesthesia

SPROTTE® Sets & IntraLong Catheter Sets




Epidural Anaesthesia

EpiLong II & EpiLong II Soft Sets




Combined Spinal Epidural Anaesthesia (CSE)

EpiSpin Sets

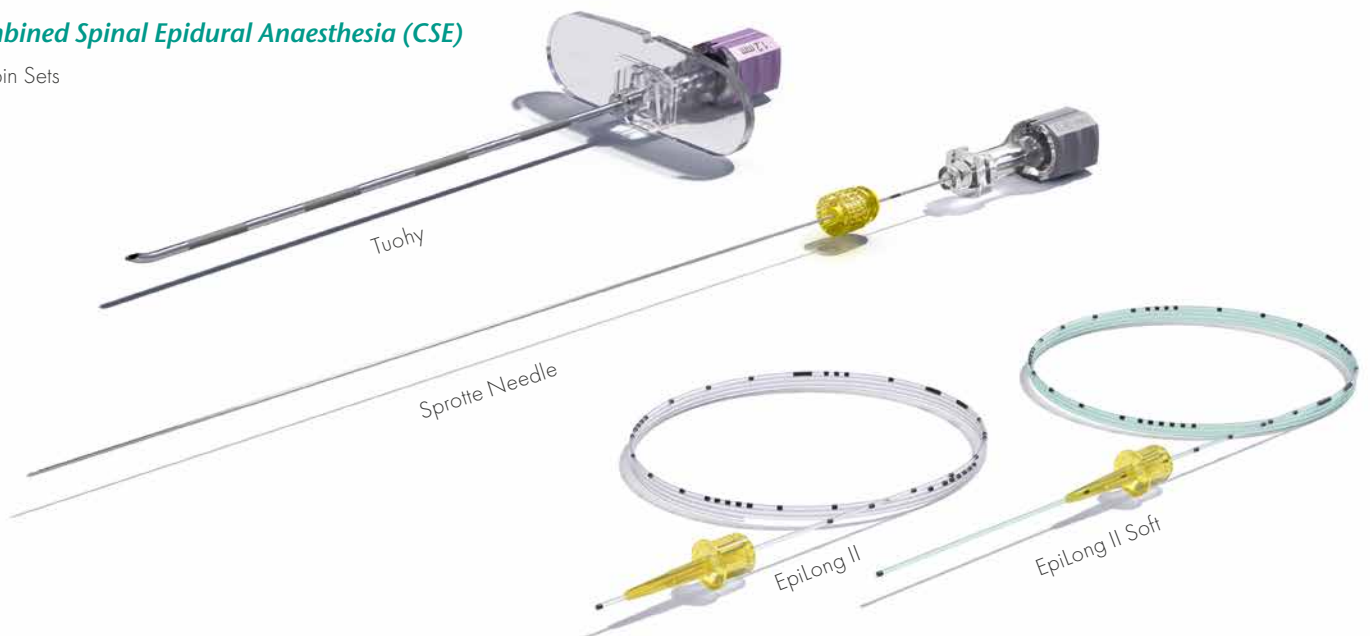


NERVE BLOCKS

Single Shot Technique

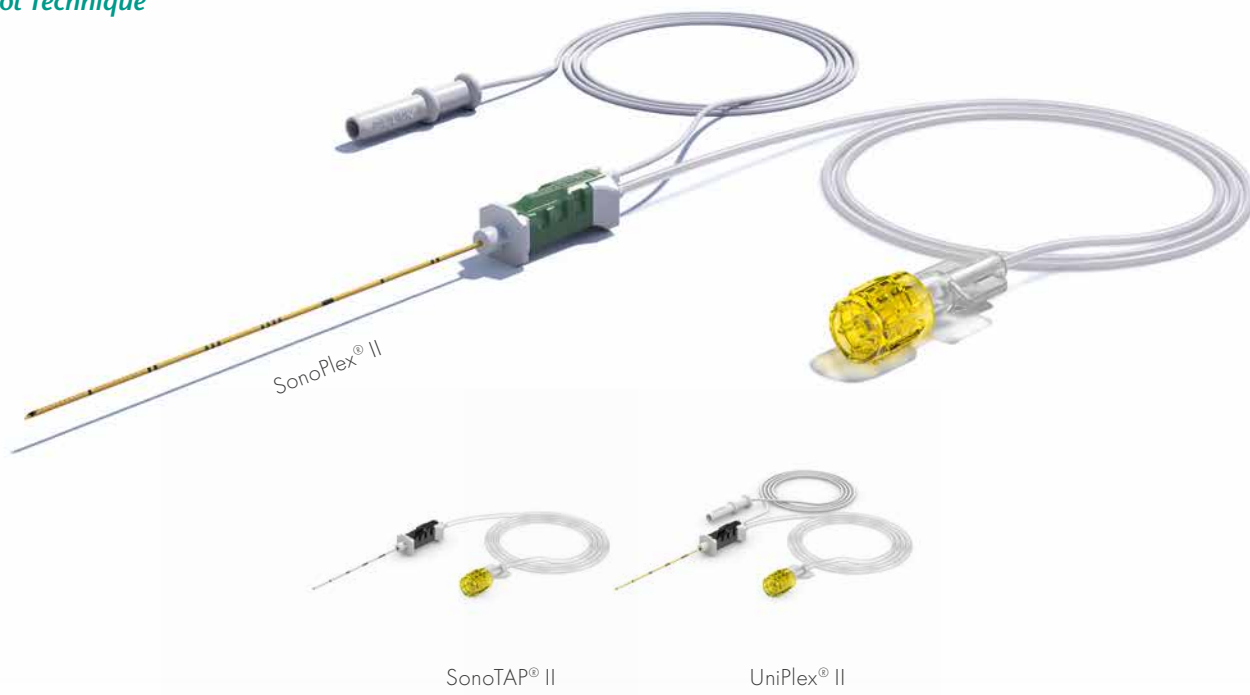

Continuous Technique

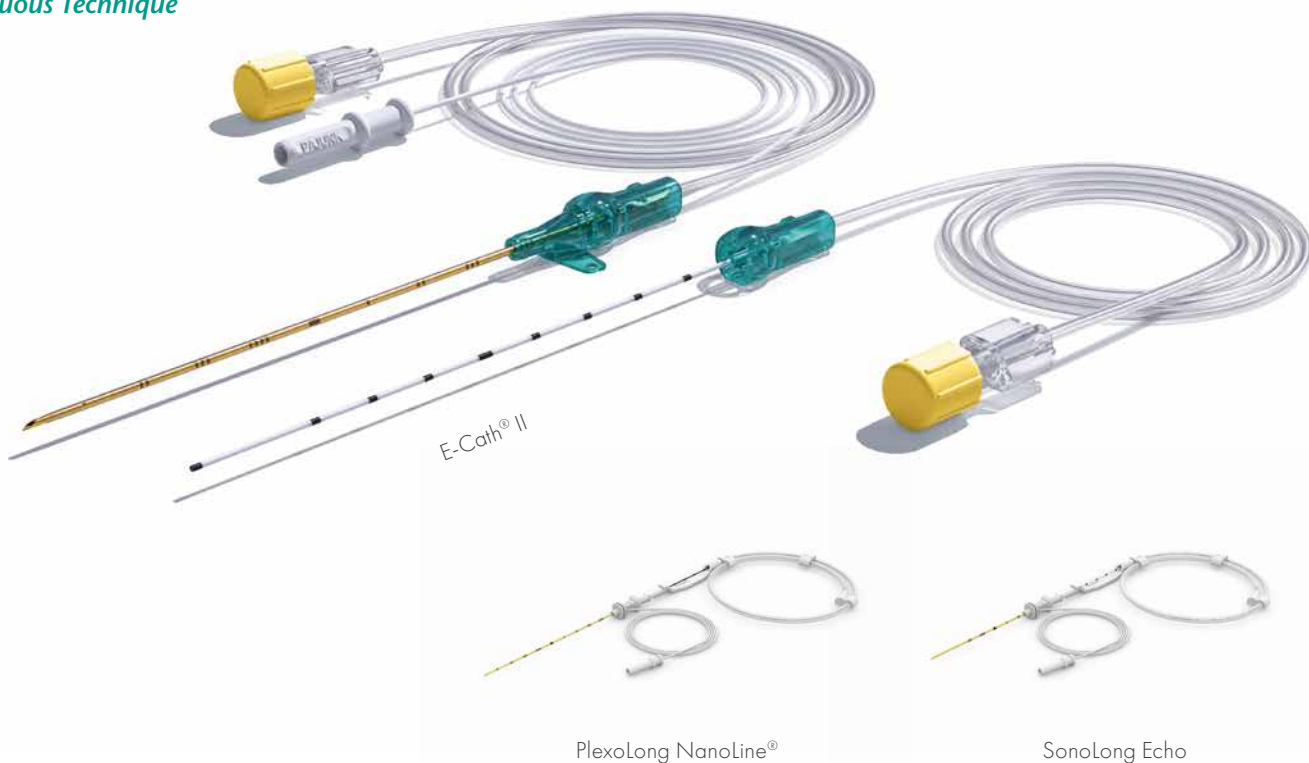

LUMBAR PUNCTURE

Lumbar Puncture



WOUND INFILTRATION

Continuous Wound Infiltration

FuserPump & InfiltraLong Sets




ACCESSORIES

Syringe Pumps Connections



Spike (vented)

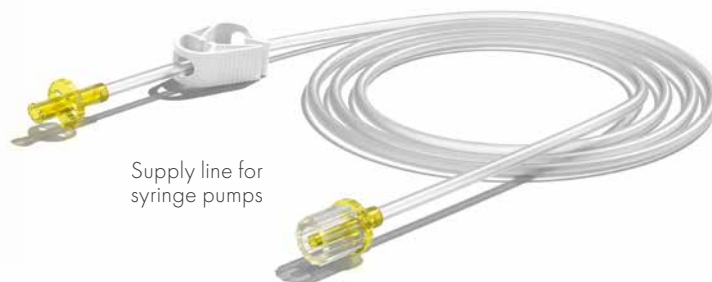

Supply line for syringe pumps

Pumps



FuserPump



Syringe 50 ml Lock



Spike (vented)

Additional NRFit Accessories



Syringe 3 ml



Syringe 5 ml



LOR 7ml Slip



Clamping Adapter



Filter Straw



Filter Needle



Draw Up Needle



FixoLong



Filter 0.2 µm



Locking Cap (male)



Locking Cap (female)

*Not all products are registered and approved for sale in all countries or regions. Indications of use may also vary by country and region. Please contact your country representative for product availability and information. Product images are for reference only.

**NRFit is a registered trademark of GEDESA used with their permission.

PAJUNK® GmbH Medizintechnologie
Global Headquarters

Karl-Hall-Straße 1
78187 Geisingen, Germany
Tel.: +49 7704 9291-0
Fax: +49 7704 9291-600
info@pajunk.com
pajunk.com

PAJUNK® Medical Produkte GmbH
Central Europe

Karl-Hall-Straße 1
78187 Geisingen, Germany
Tel.: +49 7704 8008-0
Fax: +49 7704 8008-150
info@pajunk-medipro.com
pajunk.de & pajunk.eu

PAJUNK® Swiss AG
Switzerland

Neuhofstrasse 5A
6340 Baar, Switzerland
Tel.: +41 417 609 627

info@pajunk.ch
pajunk.ch

PAJUNK® Medical System L.P.
United States of America

4575 Marconi Drive
Alpharetta, GA 30005, USA
Tel.: +1 888 972 5865
Fax: +1 678 514 3388
info@pajunk-usa.com
pajunkusa.com

PAJUNK® UK Medical Products Ltd
United Kingdom

Unit D1 The Waterfront
Goldcrest Way, Newburn Riverside
Newcastle upon Tyne, NE15 8NY, UK
Tel.: +44 191 264 7333
info@pajunk.co.uk
pajunk.co.uk